

3 IN 1 EMS Infrared air pressure therapy Machine



Warning

Dear users:

Thank you for choosing our company's delicate instrument.

In order to make full use of the product,
we sincerely advise:

- 1) Please read the user manual before any operation
- 2) Do not put any liquid onto or into the machine. Be careful about the electric connections
- 3) Avoid operations on pregnant women, or those people with epilepsy, pacemakers or any serious psychological or medical condition, or infections in the area treated
- 4) Please disassemble or assemble the equipment by professional personnel. Please check if there is any connection inside the machine dropped when some function does not work.
- 5) Clean and sterilize the spare parts after each operation.

1, Principle and Effects

This aircraft is a most popular multi-function weight-loss device, combined with functions of far-infrared fat dissolving, air-pressure lymphatic drainage, as well as low-frequency stimulation. Operate far-infrared heating after massage (by beautician), then start the lymphatic drainage operations, and at last use low-frequency to tighten skin, to achieve systematic and comprehensive effects of fat losing and size reducing.

1) Far Infrared

The far-infrared slimming treatment can help increase basic metabolic rate (BMR), which is the fastest and natural way to lose weight.

As the body warms, in addition to sweating, the blood thins and the speed at which it circulates around the body increases, while the body tries to cool itself down. The effect is to raise the temperature at the extremities of the body to the same temperature at the core of the body – a temperature at which body fat can be broken down.

2) Pressotherapy

The pressotherapy, with a set of cloth-like accessory and 20 pieces of air bags (10 pairs, or 10 pieces on each side of body), which allows for pressure by cycles or all-together for both the whole body and certain parts of the body, is a compression system designed to increase the venous and lymphatic flow and enhance extra-cellular fluid clearance.

It can help to promote the body's natural toxin clearing functions. The revitalization and oxygenation of the tissue helps to slim and redefine the legs, stomach and arms while enhancing skin tone.

It is highly recommended for the treatment of cellulite, which is a safe alternative to liposuction. It relieves pain and swelling, and provides immediate comfort. It redefines the legs, stomach and arms while enhancing skin tone.

3) EMS

By creating electronic pulse to stimulate special points on the body, EMS (Electro Muscle Stimulation) can effectively regulate bio-electric and endocrine

system and accelerate metabolism , reaching remarkable effect of weight losing and skin tightening.

2 . Installation Indication

1) Connectors on the suit



1) -1 Details



infrared connection

Please notice : press the pink part , do not loss the rubber rings.

Far infrared cables and pressotherapy tubes to be connected to machine




2) -1 details

Presso tubes to be connected to the machine , the side with Arrow is upward.



Far infrared cables to be connected to the machine.

2) Connectors on the machine

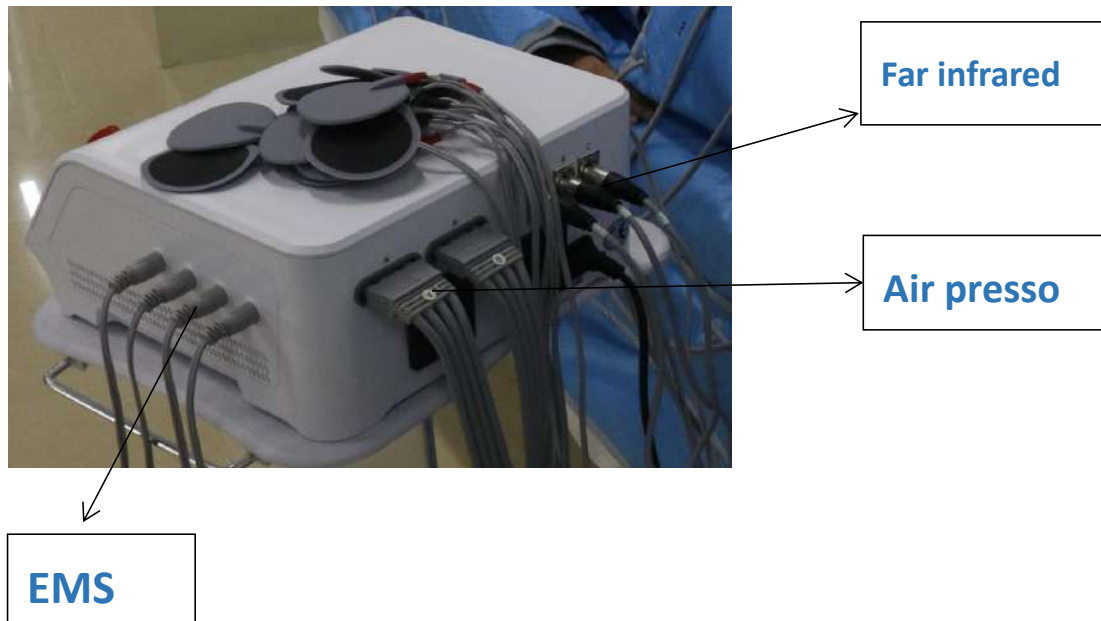

3) connection of EMS



EMS cables to be connected to machine


3 ,Operation Screen Indication

1) Connect the power supply , (please pay attention to 110V or 220Vvoltage) , turn on the machine and the default display is for pressotherapy as below



Inflation time and use for each part (the unit is second)

Pressotherapy mode indication (there are 6 modes)

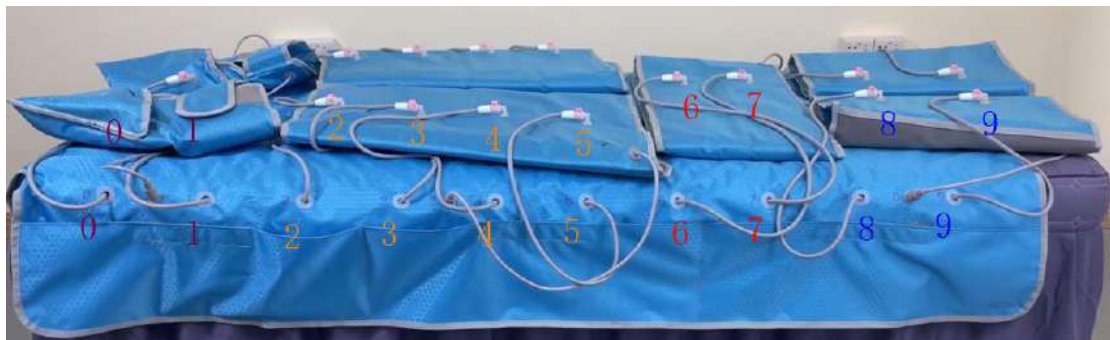

Mode 1 : it inflates from 0 to 9 . it deflates 0 when it inflates 1 , and deflates 1 when it inflates 2 , , then it inflates 0 to 9 together , and then it repeats the program .

0→1→2→3→4→5→6→7→8→9 repeat 0→1→2→3→4→5→6→7→8→9

Mode 2: working on limbs only ,it inflates from 1 to 9, (missing the body parts). it deflates 0 when it inflates 1, it deflates 2 when it inflates 3, it deflates 4 when it inflates 5, it deflates 6,7 when it inflates 8; it deflates 1 when it inflates 2, it deflates 3 when it inflates 4, it deflates 5,6,7 when it inflates 8, it deflates 8 when it inflates 9,

1→3→5→8→2→4→8→9 repeat→1→3→5→8→2→4→8→9

Mode 3: it keeps inflates 0,1,2 together, it deflates 0,1,2 when it inflates 3,4 and keeps 3,4 at the same time

0,1,2→3,4→5,6,7→8,9 repeat →0,1,2→3,4→5,6,7→8,9

Mode 4: from Limbs to body:

0,1,9→2,8→3,6,9→4,5 repeat→0,1,9→2,8→3,6,9→4,5

Mode 5: it inflates from 0 to 9, it deflates 0 when it inflates 2, and deflates 1 when it inflates 3, then it inflates 0 to 9 together, and then it repeats the program.

9→8→7→6→5→4→3→2→1 repeat→9→8→7→6→5→4→3→2→1→0

Mode 6: it inflates from 0 to 5.

2→3→4→5 →1→2→3→4→5

It is recommended to select mode 1 at the beginning, and then choose other mode according to requests

2) Press **Far Infrared** and enter the display as below

Temperature adjust for each part (the unit is DC)

Waist :42Max Celsius

Arm :39 Max Celsius

Leg :45 max Celsius



3) Press **EMS** and then enter the display as below

EMS : 5 levels max adjustable



4 , Important Notice" Please place a piece of wrap to cover the body , and then put on the suit* Please pay attention to communicate with patient about their feeling during the operation , to adjust temperature and time ." Please supply water to patients after or during the operation" Please apply some gel or cream on the body or electric pads and fix them up with Bandage.

5, Maintenance for the cloth accessory

- 1) Do not leave the cloth being heated up without patient inside
- 2) Please open the cloth to abstract heat and clean just after finishing the operation
- 3) Take care of the connections , do not pull on the pipes roughly , do not bend the pipes excessively

Specification :

Voltage :100-110v or 220-240v

Frequenc:50HZ60HZ

Power:500 W