

HOWBODY

HOWBODY H6

BODY COMPOSITION ANALYZER

EXPERIENCE THE EVOLUTION OF TECHNOLOGY

SMART · ACCURATE · CONVENIENT · EFFICIENT



www.hbec.tn

HOW?

Body —————>


Understanding Your Body = Taking Control of Yourself

What is Body Composition Analysis

The Birth of BIA Technology

Bioelectric Impedance Analysis (BIA) provides simple, intuitive, accurate, and scientific data to understand the body.

- BIA technology leverages the fact that fat tissue in the human body has low water content and poor conductivity, while muscles and internal organs have high water content and are easily conductive. By applying an alternating current to the body to measure impedance, it calculates body fat percentage and other body composition metrics.
- BIA technology is a non-invasive, safe, simple, and precise method for determining body composition, and it is currently widely applied.
- The concept of bioelectric impedance emerged at the end of the 18th century, giving it a history of over 230 years.
- Following the establishment of bioelectric theory, an increasing number of research institutions worldwide have been working on bioelectric impedance measurement technology; today, this technology is highly advanced.

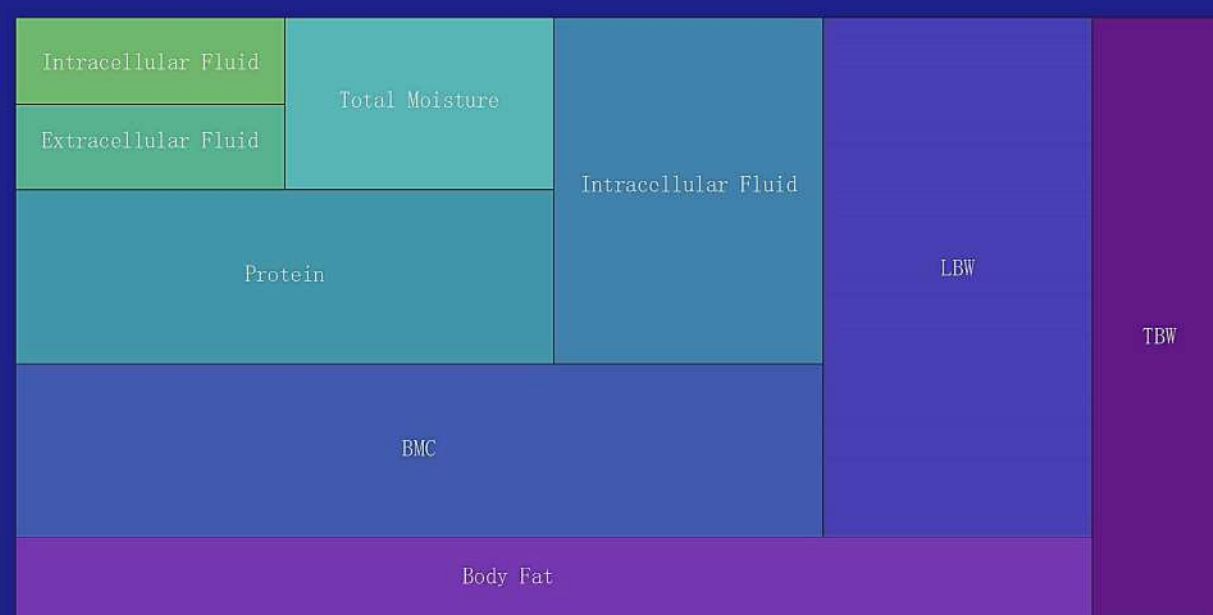


Principle of BIA Technology

(Bioelectric Impedance Analysis)



Using Bioelectrical Impedance Analysis (BIA) technology, the bioelectrical impedance is accurately measured and then applied to scientifically validated algorithms to calculate the primary components of the human body. When a safe microcurrent passes through the electrodes on the handles into the body, it easily travels through muscle, which has a high water content, but faces more resistance in fat tissue, which has higher impedance.

The closer the correlation coefficient (R^2 value is to 1, the higher the concordance between the fat mass (FM) measured by DEXA and the fat mass (FM) measured by various BIA devices. There are several methods for body composition analysis, including CT (Computed Tomography), MRI (Magnetic Resonance Imaging), DEXA (Dual-Energy X-ray Absorptiometry), and underwater weighing. Among these, DEXA is recognized as the gold standard due to its accuracy in analyzing body composition based on fat, muscle, and bone, and it is also very safe in terms of radiation exposure.

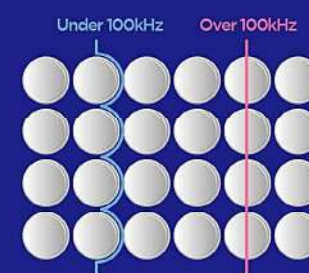

Multifrequency Measurement :

Frequencies below 100 kHz flow along the cell membrane, making them suitable for measuring extracellular water content. Fre-

ONE DEVICE SOLVES THREE MAJOR ISSUES

BODY
DATA

RISK
ASSESSMENT

CUSTOMER
SERVICE

20 BODY COMPOSITION METRICS

- HYPERTENSION
- DIABETES
- HYPERLIPIDEMIA
- FATTY LIVER

- CUSTOMER MANAGEMENT DASHBOARD
- DATA COMPARISON OVER TIME
- INCREASED CUSTOMER ENGAGEMENT




PRECISION MEASUREMENT SYSTEM

CHIP + ALGORITHM + STRUCTURE

Multifrequency Bioelectrical Impedance Analysis (DSM-BIA) Technology

The most accurate in the country (achieving a linear correlation of 0.97 with industry benchmark products)

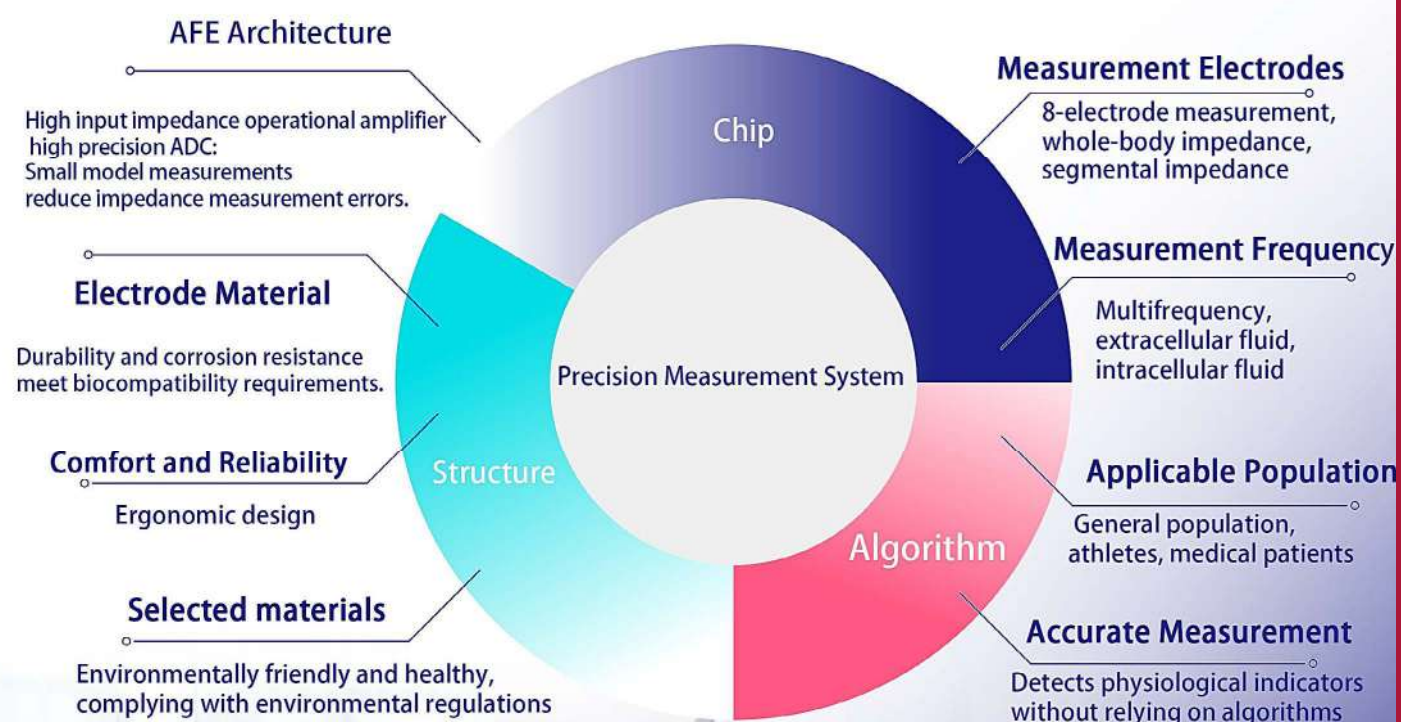

Most Accurate Technology in the Country

- Trunk resistance repeatability within 1Ω
- Limb resistance repeatability within 10Ω
- 90% within 5Ω



Breakthrough in Technical Challenges

Surpassing industry technical standards

- 8-Electrode, Multi-Segment Measurement
- Multifrequency (10K/50K/250K)
- Model construction independent of age and gender

Technical Difficulty

- Successfully Eliminated the Influence of Various Parasitic Factors
 - Significantly improved the accuracy of impedance measurement
 - Special Enhancement of Trunk Impedance Accuracy
Leading industry breakthrough
-
- The accuracy and applicability of algorithm models can extend from the general population to include professional fitness enthusiasts and individuals with hidden obesity.



LEADING ALGORITHM FOR CHILDREN




中国科学院自动化研究所
INSTITUTE OF AUTOMATION CHINESE ACADEMY OF SCIENCES



首都儿科研究所
CAPITAL INSTITUTE OF PEDIATRICS

The joint project, developed in collaboration with the China Children's Institute and the Institute of Automation, Chinese Academy of Sciences, has conducted in-depth research on measurement, data models, and other technologies for children.

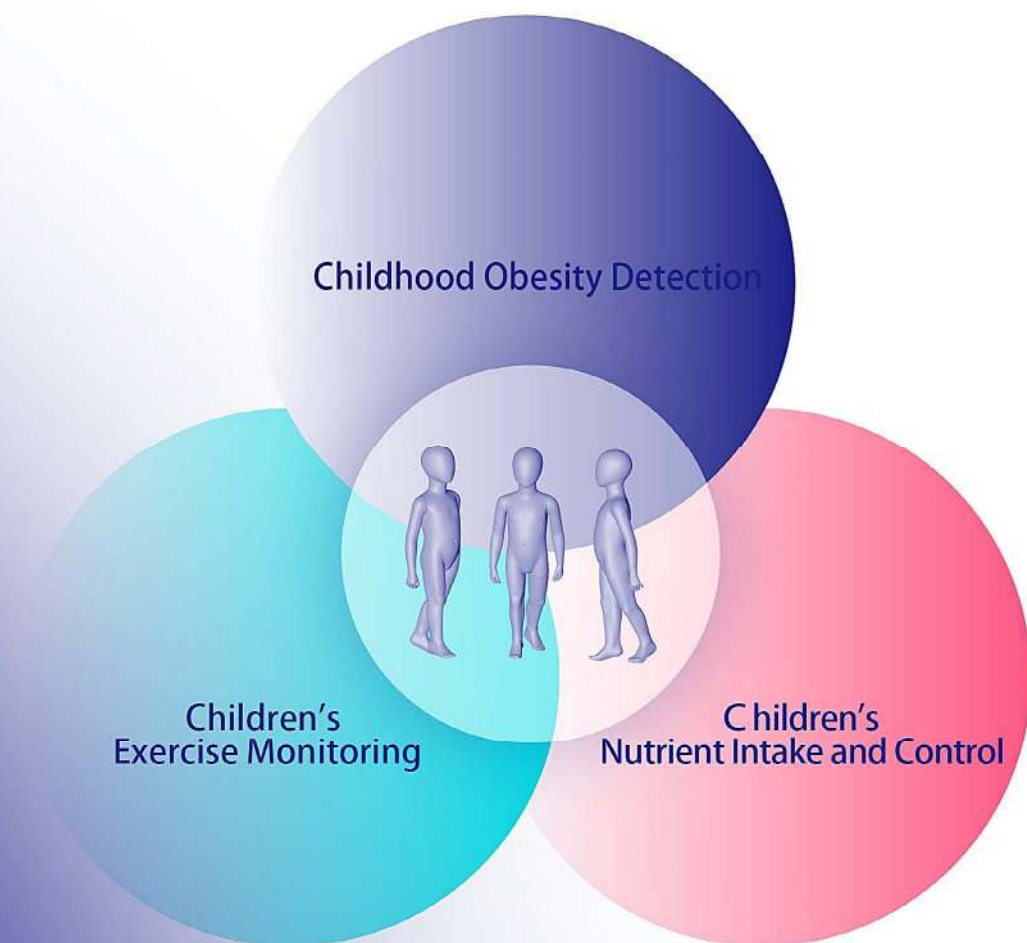

After systematic engineering optimization, the whole machine and large medical equipment (INBODY370) achieve a correlation coefficient of ≥ 0.97



Type	BMI	Body Fat Percentage (%)		
		Traditional Model	New Model	Large Medical Equipment
Fitness(Male)	24.9	20.7	16.6	15.9
Hidden(Male)	23.6	19.6	25.4	25.8
Fitness(Female)	23.3	28.9	20.3	20.5
Hidden (Female)	22.1	21.3	31.9	32.3
Children	—	14.3	12.1	12.4

Experts from the China Children's Institute and the Institute of Automation, Chinese Academy of Sciences, conducted overall planning for the project, focusing on childhood obesity detection, children's exercise monitoring, and children's nutrient intake and control.

Measurement Projects

30+

Measurement Indicators (Extensible)





Analysis Projects

10+

Measurement Indicators (Extensible)

- Human Body Composition Analysis
- Body Shape Analysis
- Obesity Analysis (BMI, Body Fat Percentage)
- Segmental Fat Analysis
- Segmental Muscle Analysis
- Muscle-Fat Analysis
- Weight Control
- Fat Control
- Muscle Control Recommendations
- Disease Risk Analysis (Hypertension, Fatty Liver, Hyperlipidemia, Diabetes)

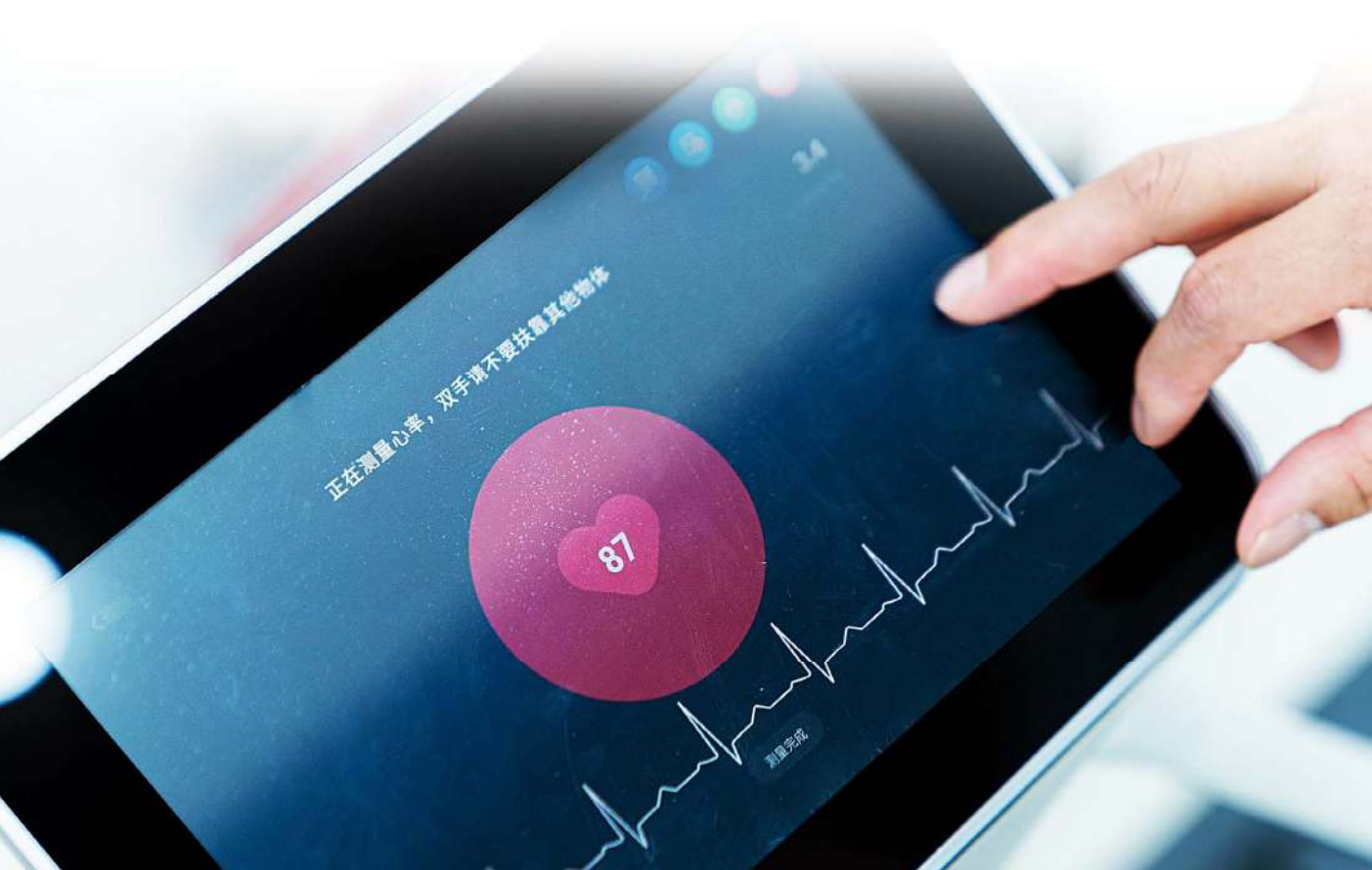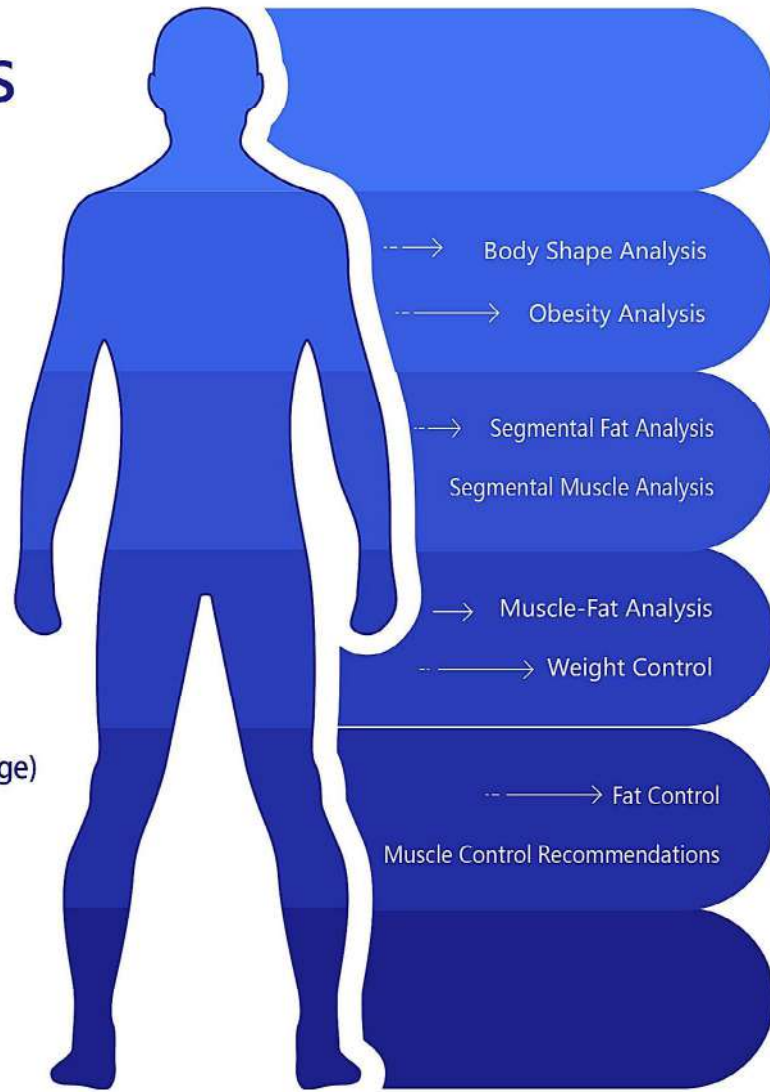

www.hbec.tn

HEALTH & BEAUTY EQUIPMENT COMPANY

HBECE

10.1-Inch
HD Screen
10-Point
Capacitive Touch

Utilizes
Android system
offering stable and
reliable performance

Multiple data
management methods
Integrated data management

Mechanical handle
with multi-angle rotation
Ensures measurement accuracy
and easy storage

Wear-Resistant
Silicone Pads
Colors can be replaced

Metal Profiles with
Skin-Like Plastic
Lightweight and simple design

Foldable from
Multiple Angles
Convenient for transportation
and maintenance

Multi-Scenario
Applications
Suitable for gyms, hospitals,
schools, beauty institutions



www.hbec.tn

HEALTH & BEAUTY EQUIPMENT COMPANY

HBECE

APPLICATION SCENARIO



GYM



HOSPITAL



SCHOOL



BEAUTY SALON



HEALTH CHECKUP




YOGA&PILATES



HEALTH RECOVERY CENTRE



CHILDREN'S PHYSICAL FITNESS

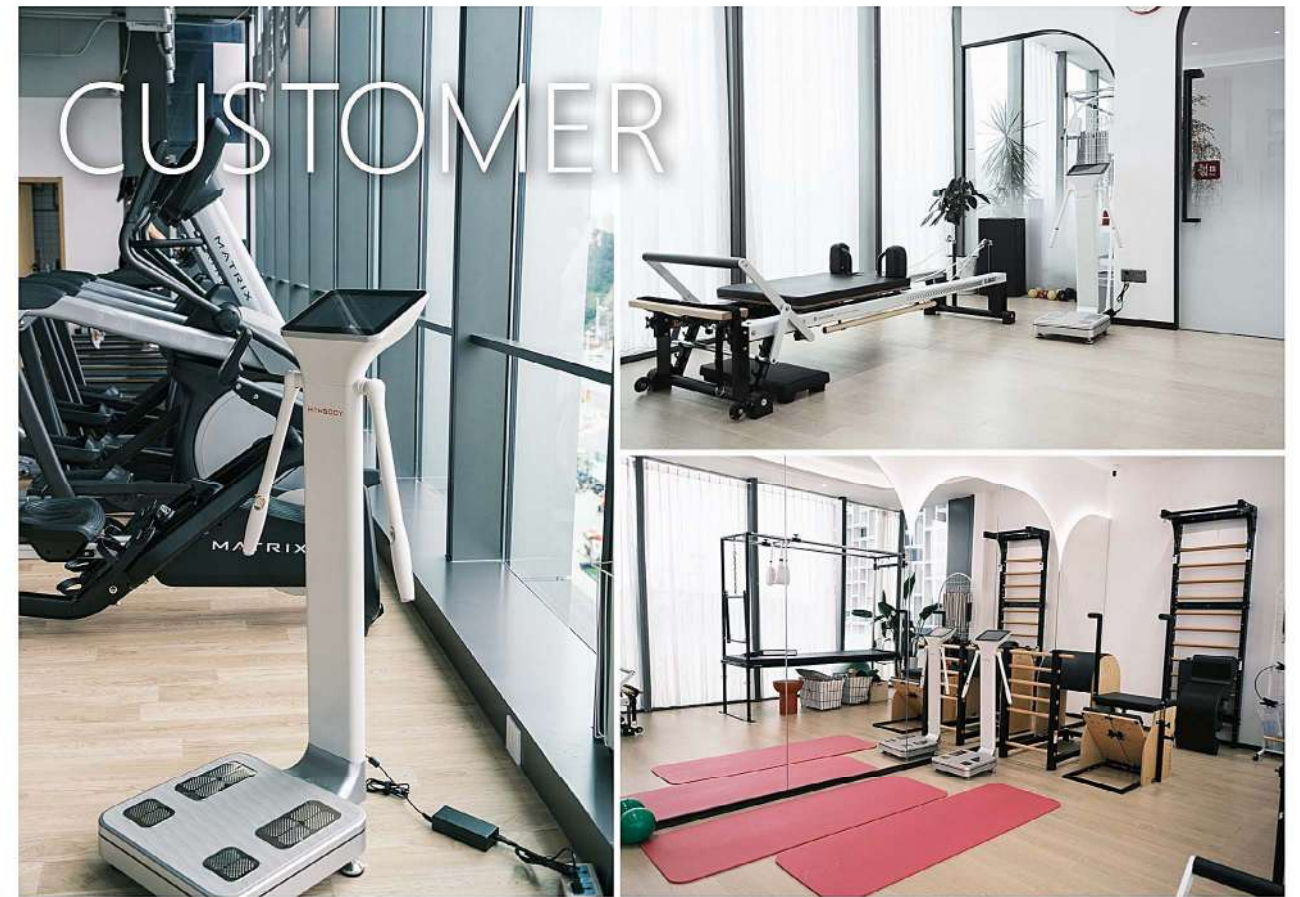

CUSTOMER




CASES

DEVICE VERSIONS



ONE DEVICE EQUALS

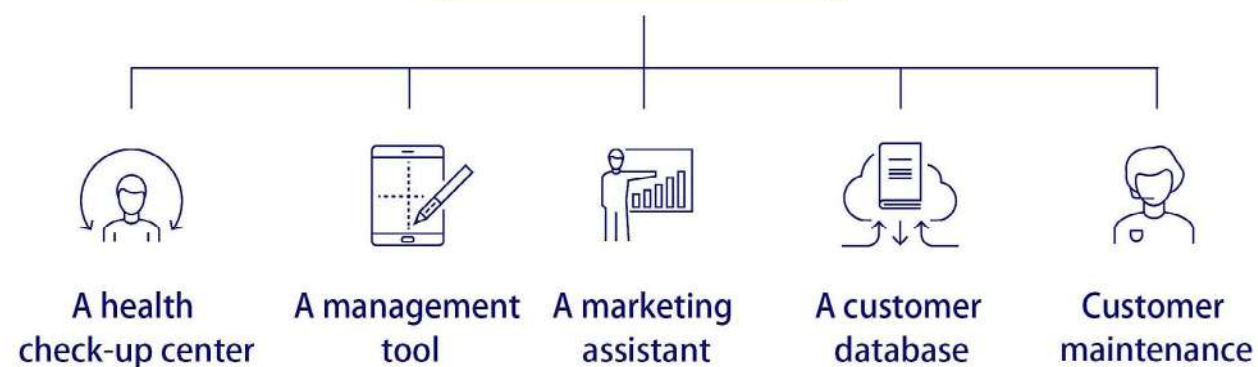

REPORT

FOUR WAYS TO VIEW THE REPORT:

Device screen



Printed report



Electronic report



Web portal

Exclusive user management platform for easy viewing and printing of historical reports.

PRODUCT SPECIFICATIONS



Measurement Electrode: 8-point contact electrode

Measurement Principle: Multi-frequency bioelectrical impedance analysis technology

Impedance Measurement: Frequencies: 10kHz, 50kHz, 250kHz

Measurement Locations: Whole body impedance、 Segmental impedance
(left arm, right arm, left leg, right leg, trunk)

Age Range for Measurement: 3 to 99 years old

Weight Range for Measurement: 5 to 200kg

Measurement Items: Protein、 Minerals、 Total body water、 Fat mass、 Muscle mass、 Fat-free mass、 Body weight、 Skeletal muscle、 BMI、 Body fat percentage、

Comprehensive score、 Basal metabolic rate、 Body shape assessment、 Heart rate

Segmental assessment (muscle mass and fat mass of left/right arms, left/right legs, trunk)

Control targets (weight/fat/muscle)、 Visceral fat level、 Waist-to-hip ratio、 Impedance values

Touch Display: 10.1-inch capacitive touch display

Power Supply: Input: AC100-240V, 50/60Hz Output: DC12V, 3A

Communication Methods: WiFi, Ethernet, Bluetooth, USB

Dimensions: 650 x 400 x 1125mm (LxWxH)

Weight: Approximately 25kg

ABOUT HOWBODY



HOWBODY specializes in detection and assessment solutions. Its primary applications are in the fields of sports and fitness, yoga and Pilates, beauty industry, education, and healthcare, with a commitment to advancing the health assessment industry.



In collaboration with a publicly listed company, HOWBODY spent three years on research and development, and in 2023, proudly launched the HowBody H5 Intelligent Body Composition Analyzer.

This device features precise measurement and data management systems, with data that aligns closely with international leading brands, achieving a linear correlation of 0.97 with industry benchmark products.

HOWBODY provides high-quality services to various sectors including fitness centers, educational institutions, and healthcare, earning unanimous praise from users.





GUANGZHOU HOWBODY HEALTH TECHNOLOGY CO., LTD



Scan

Wave Form

Cloud

Component

Model

Data

Analysis

Visualization

HEALTH & BEAUTY EQUIPMENT COMPANY

HBC

www.hbec.tn



Human Body Composition Analysis Report


ID	Measurement Time	Gender	Age	Height	Current Weight

Human Body Composition Analysis

Unit: Kilograms

Item	Measurement Value	Muscle Mass	Lean Body Weight	Weight	Standard Range
Body Water					
Protein					
Inorganic Salts					
Body Fat					

Body Type Analysis

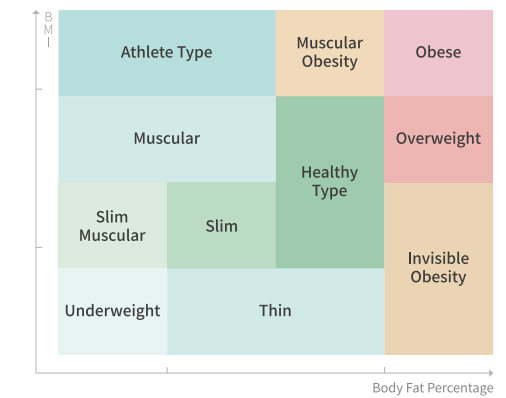


Muscle Fat Analysis

Unit: Kilograms

Item	Low Standard	Standard	Above Standard	Standard Range
Weight				
Body Fat				
Skeletal Muscle				

Segmental Fat

Unit: %

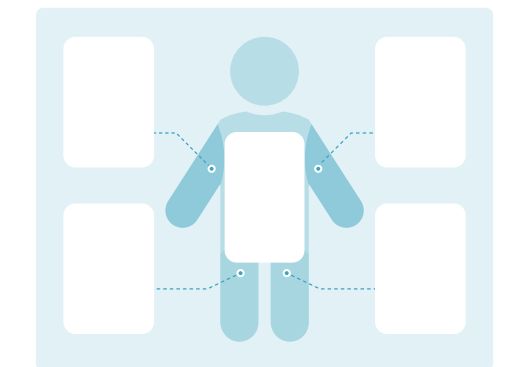


Obesity Analysis

Item	Low Standard	Standard	Above Standard	Standard Range
Body Mass Index kg/m ²				
Body Fat Percentage %				

Segmental Muscle

Unit: Kilograms

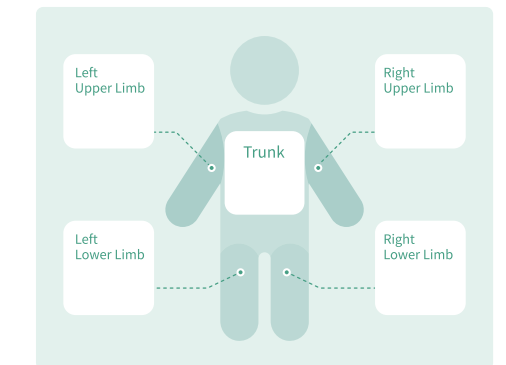


Body Data

Comprehensive Score	Basal Metabolism
Visceral Fat	Waist-Hip Ratio
Bone Mass	Body Fat Percentage

Weight Control

Weight
Muscle
Fat

Muscle Fat Analysis Standard Explanation

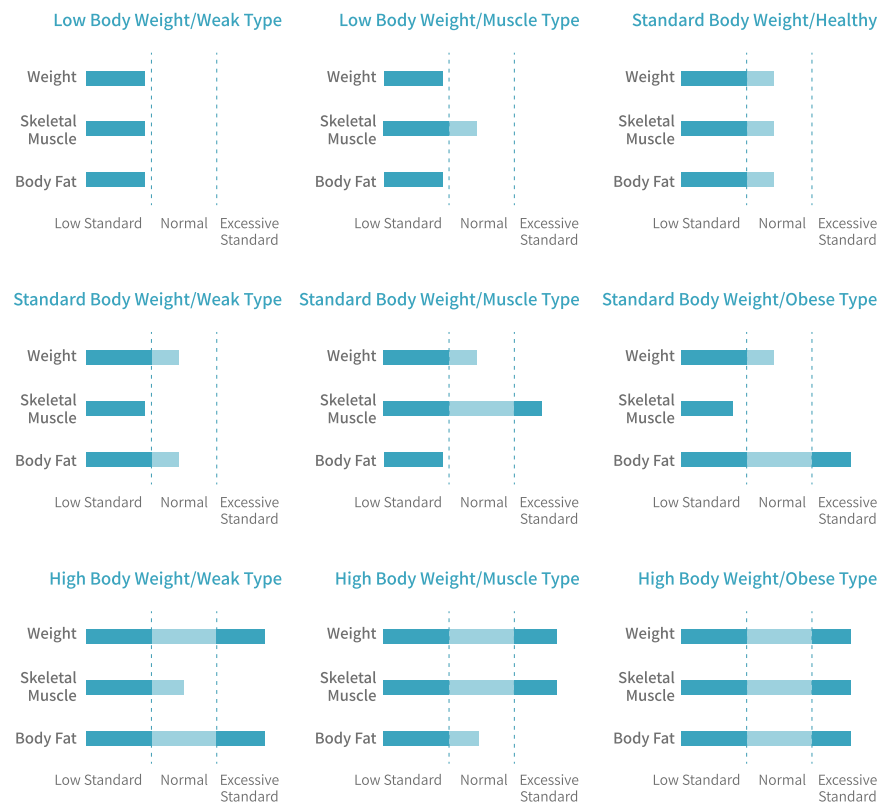

Dietary Reference

Staple Food / Seafood			Vegetables / Fruits		
Rice 166kcal/100g	drumstick 181kcal/100g	Egg 151kcal/100g	Cabbage 24kcal/100g	Chinese Cabbage 20kcal/100g	Apple 53kcal/100g
Noodles 301kcal/100g	Chicken Breast 133kcal/100g	Squid 75kcal/100g	Spinach 28kcal/100g	Cucumber 16kcal/100g	Banana 93kcal/100g
Bread 313kcal/100g	Pork 395kcal/100g	Shrimp 101kcal/100g	Cauliflower 20kcal/100g	Seaweed 13kcal/100g	Watermelon 31kcal/100g
Fish 100kcal/100g	Beef 106kcal/100g	Crab 102kcal/100g	Celery 13kcal/100g	Leeks 25kcal/100g	Orange 48kcal/100g
Beverages / Desserts			Western Food		
Cola 45kcal/100ml	Beer 32kcal/100ml		Hamburger 256kcal/100g	Steak 157kcal/100g	Sandwich 234kcal/100g
Soy Milk 31kcal/100ml	Coffee 218kcal/100g		Fried Chicken Wings 204kcal/100g	French Fries 201kcal/100g	Roast Turkey 141kcal/100g
Milk 54kcal/100ml	Milk Tea 52kcal/100ml		Pizza 235kcal/100g	Hot Dog 250kcal/100g	Salad 40kcal/100g
Ice Cream 74kcal/100g	Frappé 93kcal/100g		McNuggets 281kcal/100g	Spaghetti 351kcal/100g	Sardines 155kcal/100g

The data listed above may vary due to differences in the actual production process and added ingredients. The specific calorie values may also differ, and the data is for reference only!

Obesity Analysis Indicator Explanation

1 BMI (Body Mass Index) kg/m2

It is an indicator that determines obesity based on the ratio of weight (kg) to height (m2). Since it only considers weight and height, BMI can only assess external obesity.


2 PBF (Body Fat Percentage)


It considers not only weight but also the proportion of body fat to weight to determine obesity. Weight is not a necessary condition for determining obesity, just as having more muscle mass can lead to higher weight. In the diagnosis of obesity, body fat percentage is more important than weight.

3 WHR (Waist-Hip Ratio)

The circumference ratio between the waist and hips is usually used to determine abdominal obesity because it is closely related to it. A WHR above 0.85 for females and above 0.9 for males is considered abdominal obesity. In adults, abdominal obesity mostly leads to the accumulation of visceral fat, and excessive fat molecules increase blood concentration, leading to obesity-related diseases such as hypertension, cardiovascular diseases, and diabetes.

Explanation of Standard Range for Body Composition Analysis

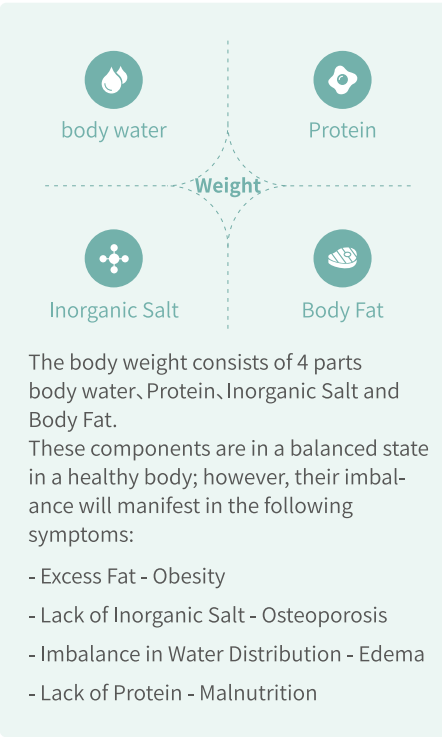

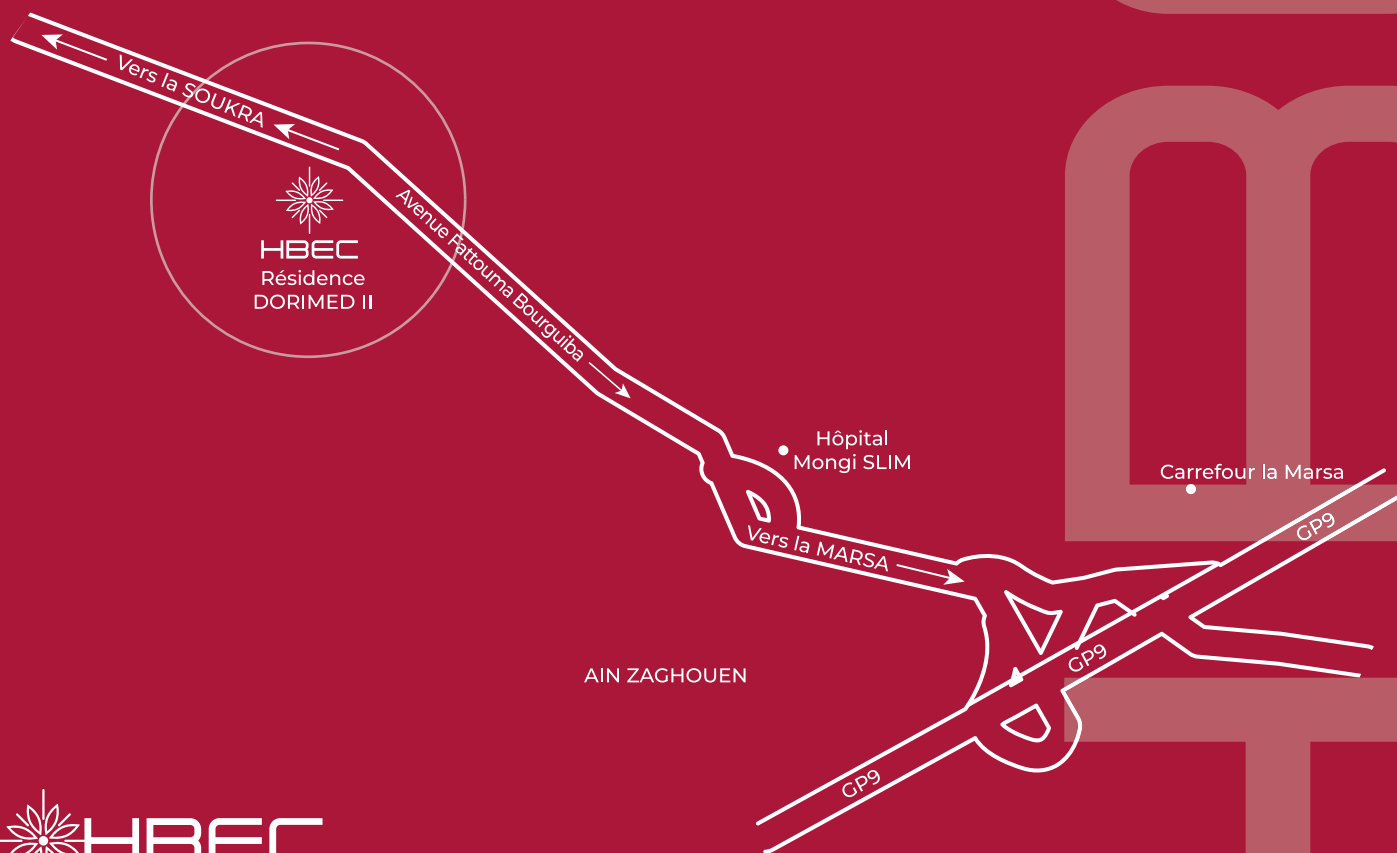

HBEC - Health Beauty Equipment Company
Bureau F 2.5 - 2^{ème} étage Résidence Dorimed 2-croisement de l'avenue Fattouma Bourguiba et de la rue d'Irak,
la Soukra 2035 - Ariana -

☎ (+216) 93 716 301 / (+216) 90 157 560

✉ contact@hbec.tn

🌐 www.hbec.tn