




7 Joint Laser Arm

7 joint laser arm from Korealab Inc.South Korea, stable output, small light loss during the output process, and uniform light energy.

USA Coherent laser

USA Coherent laser generator ensures a high level and good quality.

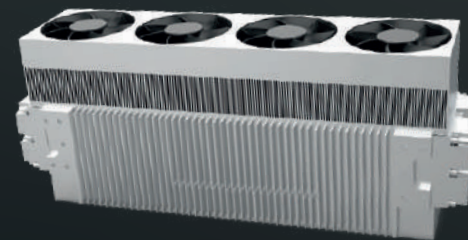

USA imported Coherent RF tube

Specification

Wavelength	10600nm
Laser Source	RF tube
Power	70W
Fractional San Shape	Square, Rectangle, Circle, Triangle, Hexagon Ellipse Line
Scan Size	20*20mm
Dwell Time	0.1ms - 10ms
Dot Pitch	0.1mm- 2.6mm
Cooling System	Air cooling
Scan Patterns	Standard, Random, Scatter
Power Supply	AC 110-220V 50/60Hz



Fractional CO2 Laser

CO2 Laser Machine

Skin Resurfacing

Acne scar removal

Freckle removal



55 157 506
90 157 560



www.gnie.com




grtarek@yahoo.fr




Imm Emmeraude de Tunis Rue Mohamed Badra
Bureau A-2-8 Montplaisir, Tunis, 1073



How Does Fractional Laser Works?

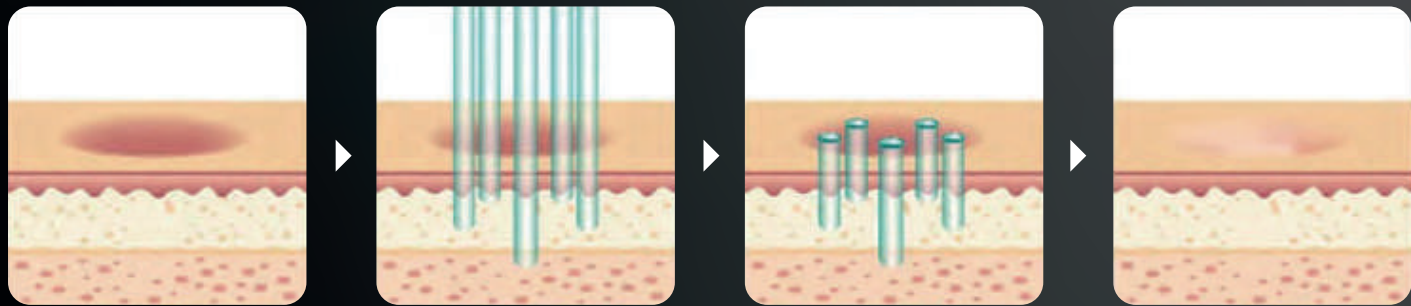


- The technology utilizes laser energy microbeams to penetrate and break down skin tissue through the epidermis.
- The laser beam will target the dermis and create tiny wounds throughout that particular area without damaging the tissues in the surrounding parts.
- Your body will naturally heal and during the process, create new, healthy tissue to replace the affected areas that have undergone treatment. Scars should be of a lighter color and your skin texture will be more even, leaving no visible wounds and allowing you to apply makeup immediately following the procedure.

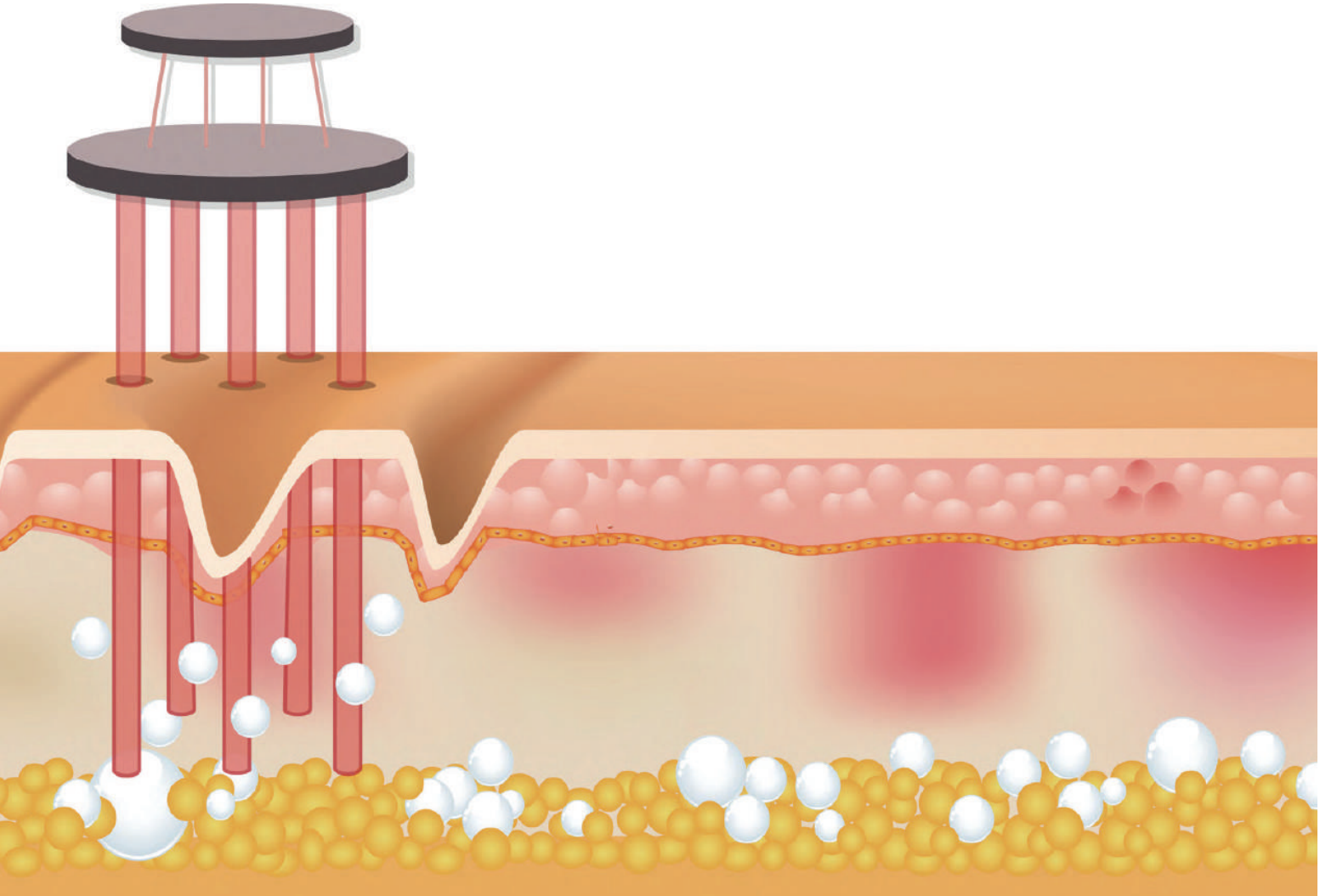

What Does CO2 Fractional Laser Do?

- ✔ Reduction of fine lines and wrinkles
- ✔ Improved skin texture
- ✔ Tightening of facial and neck skin
- ✔ Reduction of precancerous sun spots
- ✔ Improvement in the appearance of scars (e.g. acne scarring and stretch marks)

Eye bag removal

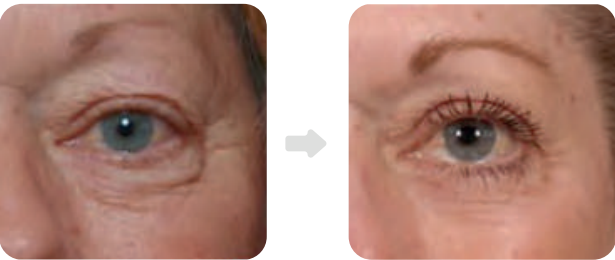


Before After

Wrinkle removal



Before After

Stretch mark removal

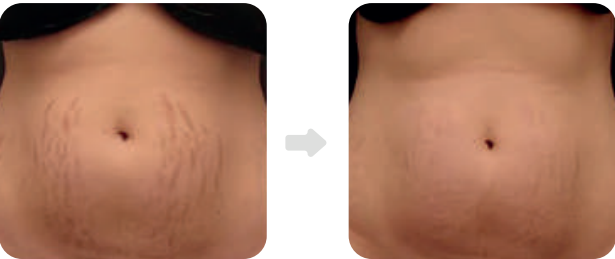

Before After

Skin resurfacing



Before After

Acne scar removal



Before After

Scar removal



Before After

Working Heads



Fractional Co2 Laser



Cutting Head



Vaginal Head

9 modes meet the treatment needs all areas



The CO2 laser's cutting head delivers focused energy to remove skin irregularities like warts and scars, offering precise tissue control and stimulating collagen for skin rejuvenation.

CO2 fractional laser causes a controlled and highly precise photothermal effect on the vaginal mucosa, promoting tissue contraction and tightening and returning its natural elasticity to the vaginal canal. The laser energy delivered along the vaginal wall heats the tissue without damaging it and stimulating the production of new collagen in the endopelvic fascia.

